

CAB-O-SIL® Fumed Silica for Pharmaceutical and Nutraceutical Applications






Cabot Corporation:

A Proud History of Innovation and Quality

Cabot Corporation has been a global manufacturer of specialty chemicals and performance products since 1882. We continue to innovate and find applications for our materials in new fields, where they help improve existing products and create exciting new offerings. Cabot is committed to the reliable delivery of pharmaceutical materials of the highest quality and consistency. Our global reach enables us to partner closely with our ever-growing list of customers in order to become their supplier of choice for quality materials, performance and service.

Current Challenges in Pharmaceutical and Nutraceutical Tablet Production

Producing tablets and capsules of uniform quality in an efficient manner can be challenging. The biggest issues tend to occur during the following stages of the production process:

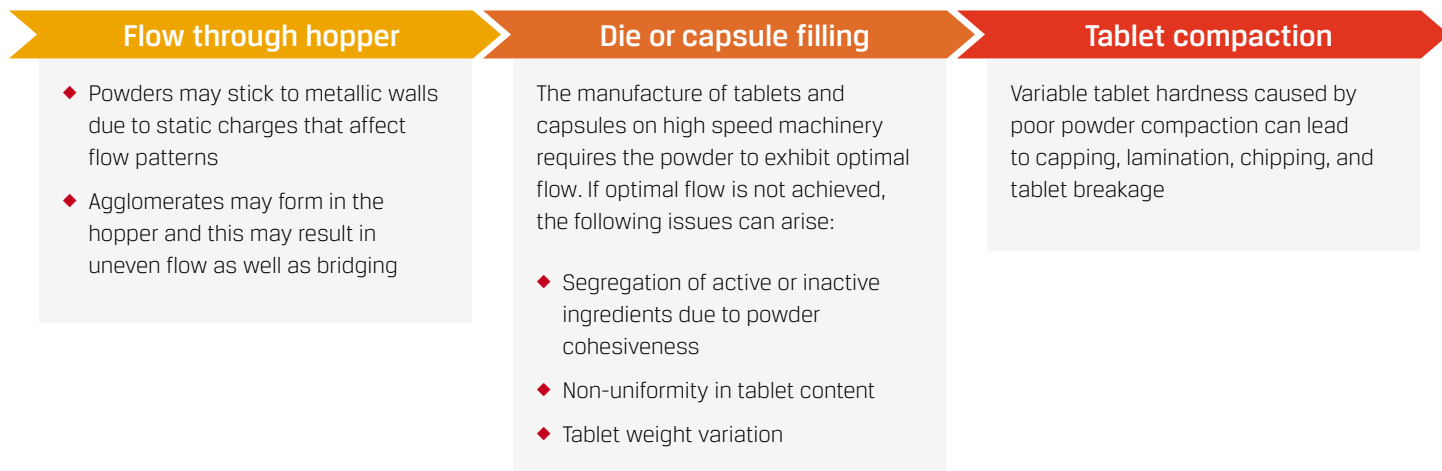

CAB-O-SIL® Product Description

CAB-O-SIL M-5P (pharmaceutical grade) fumed silica is an excipient of extremely high purity used as a multi-functional additive in the pharmaceutical and nutraceutical industries.

CAB-O-SIL fumed silica acts in the following ways to address significant issues involved in tablet production:

Flow aid	Die filling aid	Compaction aid
◆ CAB-O-SIL fumed silica's glidant and anti-static properties help improve the flow properties of powders and reduce friction and static charges in high speed tablet and capsule machines. ◆ Additionally, CAB-O-SIL fumed silica helps prevent bridging.	◆ Uniformity in tablet weight and active ingredient content result from formulations that flow evenly into the die cavity during the tableting operation.	◆ CAB-O-SIL fumed silica can aid in particle rearrangement within the die during the early stages of compression and can increase bonding strength between powder particles that form a solid tablet, helping to more effectively compact tablets/capsules. ◆ Increased tablet hardness at lower compression forces results from formulations that are effectively compacted.

Physical-Chemical Properties

CAB-O-SIL fumed silica is compatible with many pharmaceutical ingredients. It adheres readily to hydrophilic ingredients, functioning as an excellent glidant.

Pharmaceutical Product	Surface area	Bulk Density	Features
CAB-O-SIL M-5P fumed silica	200 m ² /g	<60 g/l	Easy to disperse

The data above are typical test values intended as guidance only, and are not product specifications. Product specifications are available from your Cabot representative.

Purity – CAB-O-SIL Fumed Silica versus Precipitated Silica

CAB-O-SIL fumed silica is produced through a pyrogenic process that results in one of the purest commercially available forms of synthetic amorphous silica (99.8% pure). In contrast, the production of precipitated silica requires the addition of a mineral acid to an alkaline silicate solution, resulting in a lower purity product that contains higher moisture, sulfates and chloride impurities. These impurities can impact the product in significant ways, including:

- ◆ High sulfates and chloride content can affect the stability of some active ingredients. Fumed silica has a lower level of ionic sulfates and chloride than precipitated silica.
- ◆ High levels of moisture can affect moisture-sensitive compounds. Fumed silica has a lower level of moisture than precipitated silica.

Regulatory Compliance

CAB-O-SIL® fumed silica is produced in a manner consistent with Good Manufacturing Practices appropriate for pharmaceutical excipients.


Pharmacopeia Status

Regulatory Agent	USP 37 / NF 32	EP 8 th Ed. (2014)	JP 16 th Ed. (2011)
CAB-O-SIL M-5P fumed silica	Meets all requirements of "Colloidal Silicon Dioxide"	Meets all requirements of "Colloidal Anhydrous Silica"	Meets all requirements of "Light Anhydrous Silicic Acid"

Performance

CAB-O-SIL Fumed Silica as Flow Aid

The formation of a powder bridge in a hopper of a tablet press can lead to a phenomenon known in powder flow as "bridging." The image to the right illustrates this phenomenon:

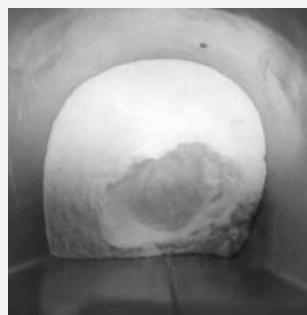

Powder Bridging in Hopper



Even Powder Flow in Hopper

How is flow measured?

The angle of repose provides a useful correlation to the flow properties of powders. The relationship between the angle of repose and powder flow is shown in Table 1 on the right.

TABLE 1: Correlation of Powder Flow to the Angle of Repose

Angle of Repose (degrees)	Flow
25 - 30	Very free flowing
30 - 38	Free flowing
38 - 45	Fair / passable flow
45 - 55	Cohesive
>55	Very cohesive

Source: Bodhmag, A. (2006). Correlation between physical properties and flowability indicators for fine powders (Unpublished doctoral dissertation). University of Saskatchewan, Saskatoon, Saskatchewan.



Figure 1 illustrates the flow properties of APAP (acetaminophen or N-acetyl-p-aminophenol), ibuprofen (isobutylphenylpropanoic acid) and aspirin (acetylsalicylic acid) tablets when tested with various concentrations of CAB-O-SIL[®] fumed silica. The results indicate that as the concentration of CAB-O-SIL M-5P fumed silica increases, the flowability of the powder mix improves, going from 'cohesive' to 'passable flow' to 'free flowing'.

FIGURE 1: Effect of CAB-O-SIL M-5P fumed silica on APAP, ibuprofen and aspirin flowability at various concentration levels

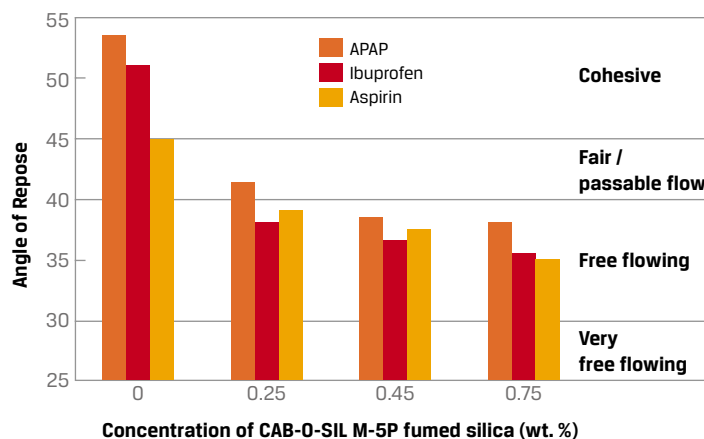


Figure 2 illustrates how the level of CAB-O-SIL M-5P fumed silica in an APAP mixture influences the flowability of the mixture. The results indicate that CAB-O-SIL M-5P fumed silica improved flow properties of the APAP mixture, taking it from a cohesive powder to fair/passable flow.

FIGURE 2: Effect of CAB-O-SIL M-5P fumed silica on APAP mixture flowability


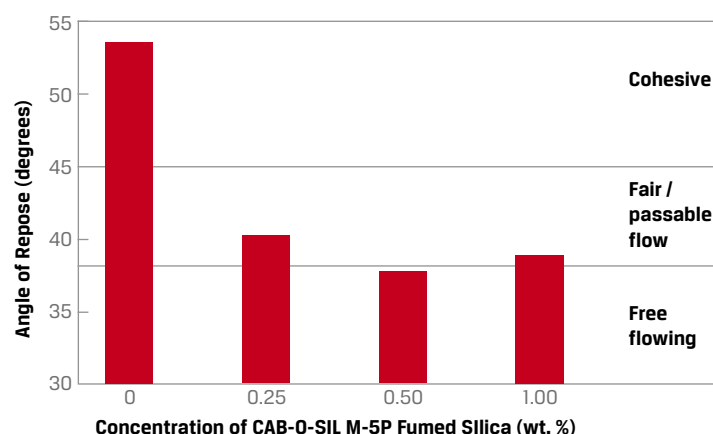


Figure 3 shows the different flow rates of Avicel[®] PH-101 Micro Crystalline Cellulose (FMC Corporation) when combined with different loading levels of CAB-O-SIL M-5P fumed silica. Note that without the addition of fumed silica, the powder does not flow.

FIGURE 3: Effect of CAB-O-SIL M-5P fumed silica on Avicel[™] PH-101 micro crystalline cellulose flow rate

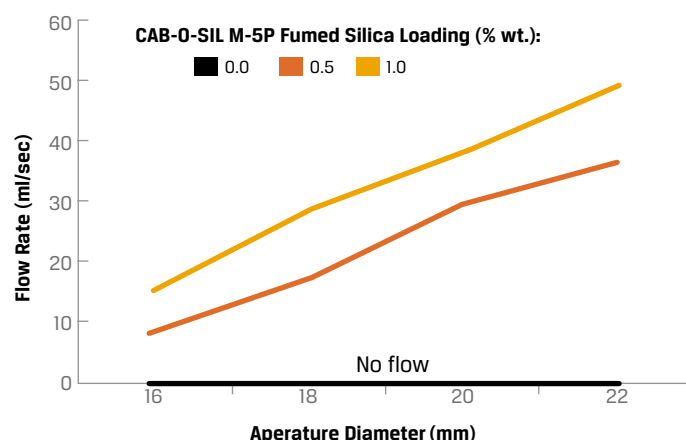

Highly Hygroscopic Materials or Materials With Static Charge Problems

For highly hygroscopic materials or materials that exhibit poor flow due to static charge, the amount of added CAB-O-SIL M-5P fumed silica may need to be 3-4 times more than that required for materials that are less prone to static.

CAB-O-SIL® Fumed Silica as a Die Filling Aid


Uniform tablet weights and uniform doses of active ingredients are dependent on the ability of the powder mix to feed into the dies in a reproducible manner. Fumed silica can be added to the formulation in order to improve the flow properties of the material and to aid particle rearrangement within the die during the early stages of die filling and compression.

Tablet Weight and Mixing Time

Formulations that flow evenly into the die cavity during the tableting operation will result in uniform tablet weight and drug content.

Figure 4 demonstrates that the weight of APAP tablets can be influenced by the concentration of CAB-O-SIL M-5P fumed silica in the formulation, as well as the time of mixing.

FIGURE 4: Effect of CAB-O-SIL M-5P fumed silica on tablet weight at various concentrations & mixing times

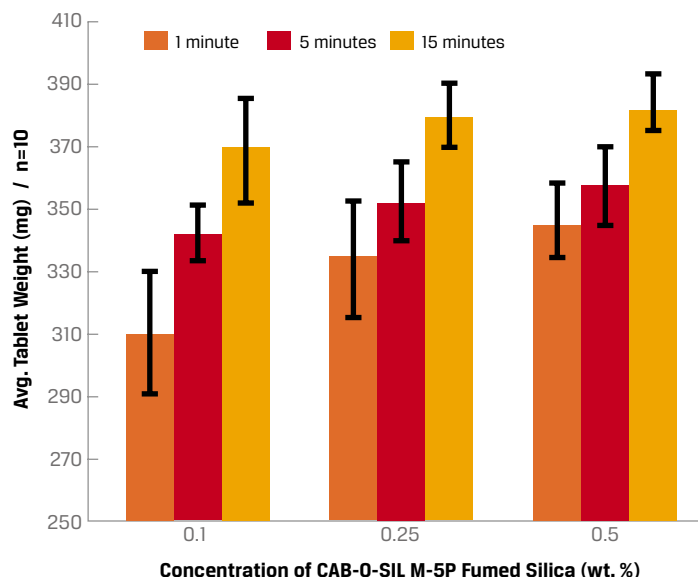


Note: The formulation includes: 10% APAP and balanced lactose. Bars correspond to minimum and maximum tablet weight.

Tablet Weight and Production Speeds

The addition of CAB-O-SIL M-5P fumed silica allows tablet manufacturers to produce high quality tablets at higher tablet press speeds than they otherwise could, therefore improving productivity.

Figure 5 illustrates that at high production speeds, the tablet weight of formulations that do not include fumed silica decreases by 11%, showcasing issues with the powder's flow and die filling. In contrast, tablets that include 0.25% CAB-O-SIL fumed silica preserve a more consistent weight under higher press speeds.


FIGURE 5: Effect of CAB-O-SIL M-5P fumed silica & press speed on tablet weight

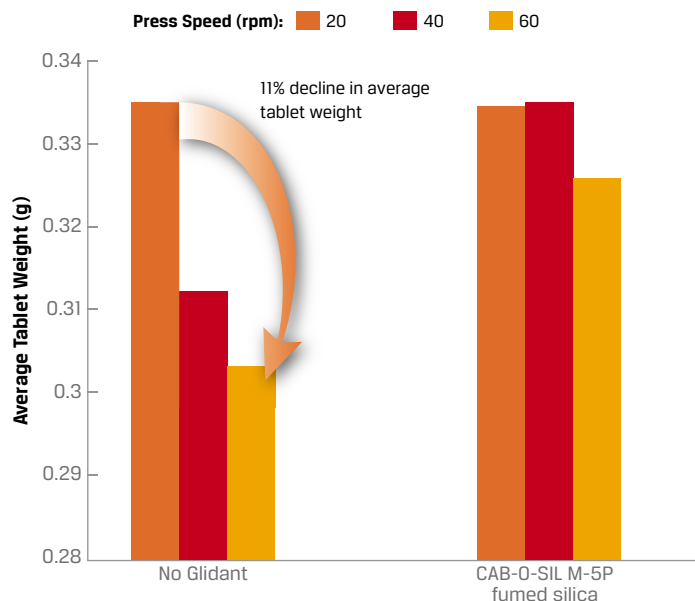

Note: Ibuprofen tablets with 0.25% CAB-O-SIL fumed silica.

Figure 6 illustrates that at high production speeds, the addition of 0.25% CAB-O-SIL® M-5P fumed silica reduced tablet weight variation anywhere from 30% to 40% when compared to the same powder mix without a glidant.

FIGURE 6: Effect of CAB-O-SIL M-5P fumed silica on ibuprofen tablet weight variation

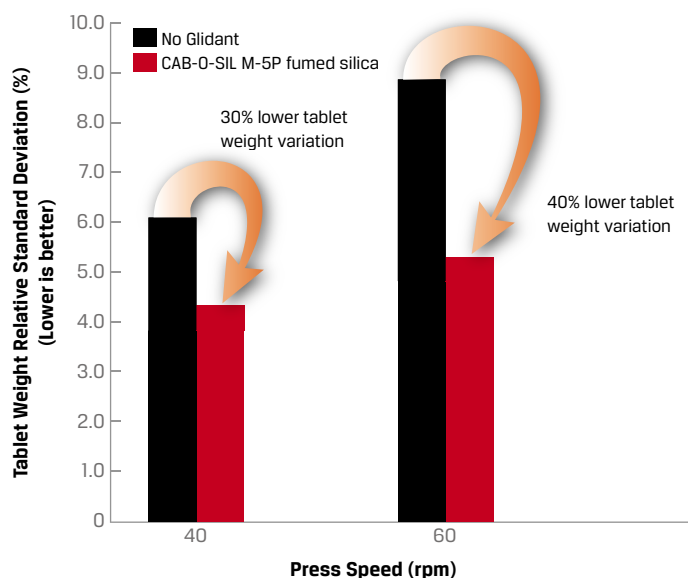


Note: Ibuprofen tablets with 0.25% CAB-O-SIL fumed silica.

API Content Uniformity

Addition of fumed silica reduces segregation of active and inactive ingredients, promoting content uniformity.

The data in **Figure 7** demonstrates that the presence of CAB-O-SIL M-5P fumed silica in a tablet formulation can decrease the coefficient of variation of active pharmaceutical ingredient (API) content in a batch of tablets. Very low variations in drug content were observed with APAP tablet formulations containing 0.5% CAB-O-SIL M-5P fumed silica after blending for 15 minutes. Tablets without fumed silica mixed for 15 minutes exhibited a coefficient of variation of 6%, which is 12 times higher than the variation seen in tablets that did contain CAB-O-SIL fumed silica.

FIGURE 7: Effect of CAB-O-SIL M-5P fumed silica on coefficient of variation of API in tablet formulations

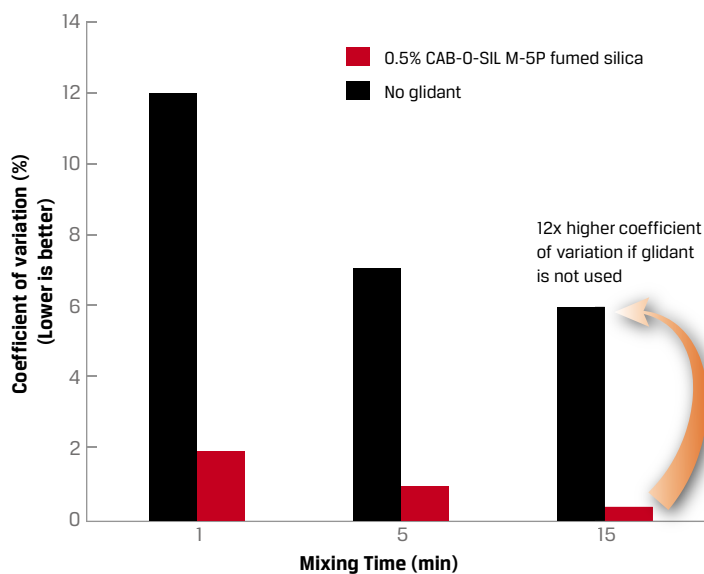

Note: These tablets contained 10% APAP, 89.5% lactose and 0.5% CAB-O-SIL M-5P fumed silica. The tablet formulation without fumed silica contained 90% lactose and 10% APAP.

CAB-O-SIL® Fumed Silica as Compaction Aid


Tablet Hardness

CAB-O-SIL fumed silica was added to ibuprofen tablets in order to test its influence on tablet hardness. The results of this testing under various compression force levels indicate the following conclusions and are illustrated in **Figure 8** below:

- ◆ The addition of only 0.25% of CAB-O-SIL fumed silica at 0.25mt of compression force facilitates an 81% improvement in tablet hardness compared to tablets formulated without a glidant. Adding 0.25% CAB-O-SIL fumed silica increases tablet hardness by 104% when compression forces reach 0.5mt.
- ◆ The addition of CAB-O-SIL fumed silica can double tablet hardness when compared to tablets formulated without a glidant.

Figure 8 indicates that the highest tablet hardness is obtained at 1.5 metric ton of compression. This is relevant, because as the compression forces are reduced, the wear on press and dyes is reduced as well.

FIGURE 8: Effect of CAB-O-SIL fumed silica on tablet hardness at various compression force levels

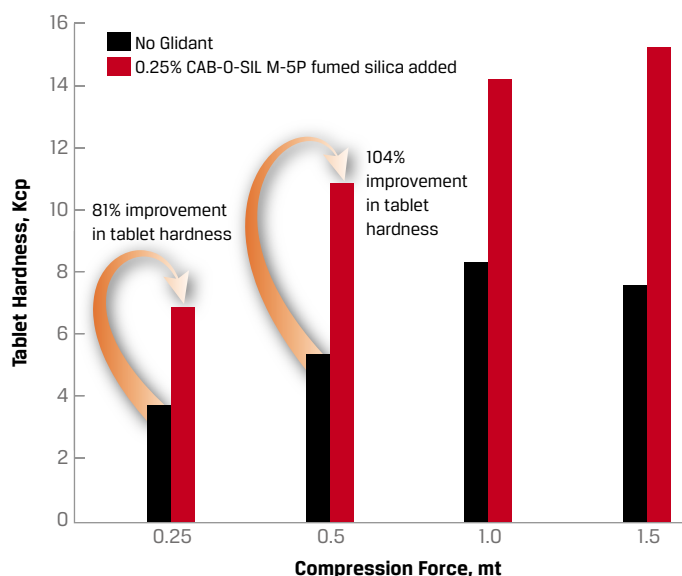

NORTH AMERICA

Cabot Corporation Business and Technical Center
157 Concord Road
Billerica, MA 01821
USA
Tel: +978 663 3455
Fax: +978 663 5471
Tel: +678-297-1300 (Customer service)

SOUTH AMERICA

Cabot Brasil Industria e Comercio Ltda.
Rua do Paraíso 148 - 5 andar
04103-000 São Paulo, SP Brazil
Tel: +55 11 2144 6400
Fax: +55 11 3253 0051

EUROPE

SIA Cabot Latvia
101 Mukusalas Street
Riga, LV-1004, Latvia
Tel: +371 6705 0984
Tel: +371 6705 0900
Fax: +371 6705 0985

MIDDLE EAST & AFRICA

Cabot Dubai
P.O. Box 17894
Jebel Ali Free Zone
LOB 15, Office 424
Dubai
United Arab Emirates
Tel: +971-4-8871800
Fax: +971-4-8871801

ASIA PACIFIC

Cabot China Ltd.
558 Shuangbai Road
Shanghai 201108
China
Tel: +86 21 5175 8800
Fax: +86 21 6434 5532

JAPAN

Cabot Specialty Chemicals Inc.
Sumitomo Shiba-Daimon Bldg. 3F
2-5-5 Shiba Daimon, Minato-ku
Tokyo 105-0012
Japan
Tel: +81 3 6820 0255
Fax: +81 3 5425 4500

The CAB-O-SIL name is a registered trademark of Cabot Corporation. The AVICEL name is a registered trademark of FMC Corporation. The FLODEX name is a trademark of Hanson Research Corporation.

The data and conclusions contained herein are based on work believed to be reliable; however, Cabot cannot and does not guarantee that similar results and/or conclusions will be obtained by others. This information is provided as a convenience and for informational purposes only. No guarantee or warranty as to this information, or any product to which it relates, is given or implied. CABOT DISCLAIMS ALL WARRANTIES EXPRESS OR IMPLIED, INCLUDING MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE AS TO (i) SUCH INFORMATION, (ii) ANY PRODUCT OR (iii) INTELLECTUAL PROPERTY INFRINGEMENT. In no event is Cabot responsible for, and Cabot does not accept and hereby disclaims liability for, any damages whatsoever in connection with the use of or reliance on this information or any product to which it relates.